

Nanoscale

Accepted Manuscript

This article can be cited before page numbers have been issued, to do this please use: Y. H. Ou, S. Zou, W. J. Goh, S. Y. Chong, G. Venkatesan, M. G. Wacker, G. Storm, J. Wang, B. Czarny, G. Pastorin and E. C. Y. Woon, *Nanoscale*, 2020, DOI: 10.1039/D0NR03784B.



This is an Accepted Manuscript, which has been through the Royal Society of Chemistry peer review process and has been accepted for publication.

Accepted Manuscripts are published online shortly after acceptance, before technical editing, formatting and proof reading. Using this free service, authors can make their results available to the community, in citable form, before we publish the edited article. We will replace this Accepted Manuscript with the edited and formatted Advance Article as soon as it is available.

You can find more information about Accepted Manuscripts in the [Information for Authors](#).

Please note that technical editing may introduce minor changes to the text and/or graphics, which may alter content. The journal's standard [Terms & Conditions](#) and the [Ethical guidelines](#) still apply. In no event shall the Royal Society of Chemistry be held responsible for any errors or omissions in this Accepted Manuscript or any consequences arising from the use of any information it contains.

Communication

Micro Cell Vesicle Technology (mCVT): A Novel Hybrid System of Gene Delivery for Hard-To-Transfect (HTT) Cells

Received 00th January 20xx,
Accepted 00th January 20xx

Yi-Hsuan Ou^a, Shui Zou^a, Wei Jiang Goh^{a,b}, Suet Yen Chong^{c,d}, Gopalakrishnan Venkatesan^a, Matthias G. Wacker^a, Gerrit Storm^c, Jiong-Wei Wang^{*c,d,e}, Bertrand Czarny^{*f}, Giorgia Pastorin^{*a,b}, Esther C. Y. Woon^{*a,g}

DOI: 10.1039/x0xx00000x

A hybrid gene delivery platform, micro Cell Vesicle Technology (mCVT), produced from the fusion of plasma membranes and cationic lipids, is presently used to improve the transfection efficiency of Hard-To-Transfect (HTT) cells. The plasma membrane components of mCVTs impart specificity in cellular uptake and reduce cytotoxicity in the transfection process, while the cationic lipids complex with the genetic material and provide structural integrity to mCVTs.

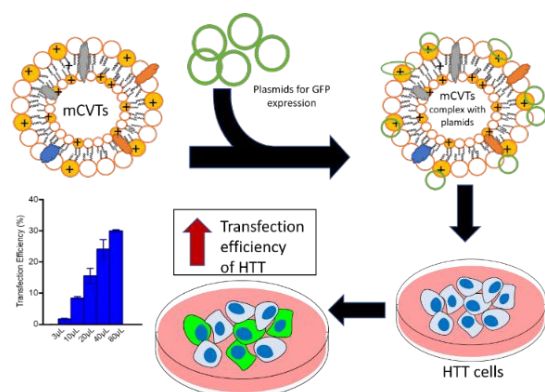


Figure 1. Schematic overview of transfection of HTT cells using mCVTs.

^a Department of Pharmacy, National University of Singapore, Singapore

^b NUS Graduate School for Integrative Sciences and Engineering, Centre for Life Sciences (CeLS)

^c Department of Surgery, Yong Loo Lin School of Medicine, National University of Singapore, Singapore

^d Cardiovascular Research Institute (CVRI), National University Heart Centre Singapore (NUHCS) and National University Health System (NUHS)

^e Department of Physiology, Yong Loo Lin School of Medicine, National University of Singapore, Singapore

^f School of Materials, Science and Engineering & Lee Kong Chian School of Medicine (LKMedicine), Nanyang Technological University, 50 Nanyang Avenue, Singapore 639798, Singapore

^g School of Pharmacy, University College London, 29-39 Brunswick Square, London WC1N 1AX, United Kingdom

*Co-first Authors, Authors contributed equally to this paper
Correspondence:

*Dr Jiong Wei Wang, National University of Singapore, Department of Surgery, Yong Loo Lin School of Medicine, 1E Kent Ridge Road, NUHS Tower Block, Level 8, Singapore 119228; Email: surwang@nus.edu.sg

*Dr Bertrand S. Czarny, Nanyang Technological University, School of Materials, Science and Engineering & Lee Kong Chian School of Medicine (LKMedicine), 50 Nanyang Avenue, Singapore 639798; Email: bczarny@ntu.edu.sg

*A/Prof Giorgia Pastorin, National University of Singapore, Lower Kent Ridge Road, 3 Science Drive 2, Singapore 117543; Email: phapg@nus.edu.sg

*Dr Esther C. Y. Woon, National University of Singapore, 18 Science Drive 4, Singapore 117543; Email: esther.woon@nus.edu.sg; School of Pharmacy, University College London, 29-39 Brunswick Square, London WC1N 1AX, United Kingdom
Electronic Supplementary Information (ESI) available: [details of any supplementary information available should be included here]. See DOI: 10.1039/x0xx00000x

Introduction

Gene delivery is the process of introducing foreign DNA into eukaryotic cells, such as the delivery of therapeutic or functional genes for alterations of cellular processes in both basic research (e.g. in drug discovery and in mechanistic studies of functional proteins¹) and clinical applications (e.g. gene therapy). Effective gene transfection requires not only the successful introduction of foreign nucleic acids into the cells but also their subsequent trafficking in and expression by the cells. This is a challenging task as DNA molecules are inherently anionic, thus poorly internalized by cells in their native form.² There is therefore tremendous interest in the development of effective gene delivery systems. Over the last decades, a great variety of carrier systems has been explored, ranging from virus-based delivery systems to non-viral vehicles.^{2–7} Despite their relatively low transfection activity, non-viral gene delivery vehicles are expected to be preferred in terms of safety and clinical acceptance, with the cationic lipid systems⁸ and polymer-based platforms^{4–6,9,10} being the prime examples. **Error! Reference source not found.** summarizes the main advantages and disadvantages of currently available transfection reagents/methods.^{1,10} An ideal non-viral gene delivery system should be able to efficiently deliver the required DNA molecules into specific cell types with minimal cytotoxicity to ensure maximal expression of the gene of interest¹¹. However, to the best of our knowledge, current commercially-available transfection reagents fall considerably short of these criteria, especially when they are used for transfect Hard-To-Transfect (HTT) cells.^{12,13}

HTT cells typically include primary cells, embryonic stem cells, and some immortalized cell lines,¹⁴ such as 3T3-L1 (murine preadipocyte),¹⁵ U937 (human monocyte),^{16,17} Jurkat (human T lymphocyte),^{6,7,9,18} HUVEC (human umbilical vein endothelial cells) and many others.^{13,19–22} **Error! Reference source not found.** shows the transfection efficiency of HTT cells with plasmid (as an example) using commercially-available reagents. One of the reasons why HTT cells demonstrate poor uptake of foreign substances and plasmid relates to their low metabolic activity.²¹ Furthermore, these cells are usually highly sensitive to cellular perturbations evoked by the transfection reagents used and potentially leading to cytotoxicity and, eventually, cell death.^{8,23}

COMMUNICATION

Journal Name

Recently, cell-derived delivery systems have gained considerable attention, in view of their potential to deliver active cargo molecules to target sites with some extent of specificity.²⁴ Indeed, a number of studies have explored the use of monocyte-derived drug delivery systems (DDS) in targeted delivery of chemotherapeutics^{25,26} and plasmid DNA (pDNA)²⁷ to tumors, or erythrocyte-derived DDS for prolonged blood circulation and targeted delivery of pDNA to the blood cells.²⁸ It is postulated that the preservation of the cellular membrane components from the parent cells would confer to these cell-based DDS the ability to recognize specific target cells, possibly promoted by the ability to evade recognition by the body's immune system.²⁹

Despite significant therapeutic interest, to date, the full potential of cell-based DDS in gene delivery remains largely unexplored. The key challenges to address include tedious and low-throughput isolation/production processes of cell-derived DDS (e.g. in the form of exosomes) and suboptimal gene loading/complexation techniques^{30,31}. Thus, while these systems seem advantageous in terms of intrinsic targeting properties, they have strong limitations regarding production yields and pDNA loading efficiencies. Hence, in order to overcome these limitations, we propose here the development of a novel hybrid gene delivery system that consists of both cellular membrane components for enhanced internalization by target cells and synthetic lipid components for enhanced pDNA loading.

In this study, we describe for the first time the development of micro-Cell Vesicle Technologies (mCVTs), a novel hybrid delivery platform obtained through the fusion of cationic lipids with cellular membranes (

Figure 2. Schematic of mCVTs and their structural components.

Figure 2). To produce mCVTs, Cell Ghosts (CG) are first obtained by treating cells with a hypotonic aqueous solution, as described by Goh et al.³¹ The hypotonic solution creates transient openings in the plasma membranes to allow the removal of the intracellular contents of the cells, while maintaining the integrity of the plasma membranes. CG derived from different cells suspended in aqueous medium are used to hydrate a dry thin film of DOTAP lipids and subsequently extruded through polycarbonate filters to provide mCVTs with a similar mean hydrodynamic size and zeta potential as compared to extruded DOTAP liposomes. In this study, we demonstrate that mCVTs are able to complex with plasmids and be internalized by target cells with minimal cytotoxicity. Interestingly, when mCVTs produced from different cells were compared *in vitro*, they showed significantly higher transfection efficiency for cell types from which the mCVTs were derived from.

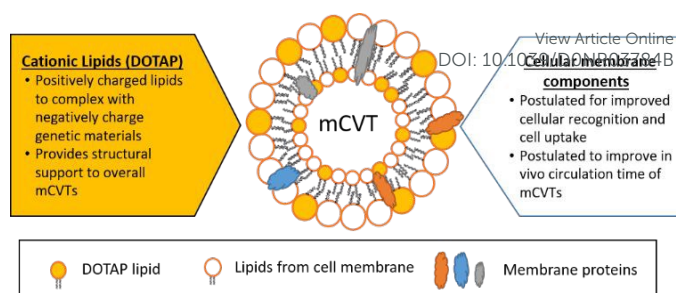


Figure 2. Schematic of mCVTs and their structural components.

Results & Discussion

Production of Cell Ghosts (CG) to ensure reproducibility and safety profile of mCVTs

Removal of intracellular contents and lack of viability are important aspects for the safe use of CG and for reducing any incompatibility and interference in subsequent applications. Cell emptying plays a vital role in mCVTs production, as cells that are not sufficiently emptied have greater aggregation propensity that may result in variation in mCVTs' size and eventually affect the reproducibility among different batches.

CG were produced by varying the osmotic pressure of the medium in which the cells were suspended (**Figure 3A**). Harvested 3T3-L1 cells (fibroblast-like mouse preadipocyte cells) were incubated in PBS hypotonic solution with increased osmotic pressure to facilitate the formation of transient pores in the cell plasma membrane; at these conditions, the cytoplasm and organelles (e.g. nuclei, ribosomes, mitochondria) were postulated to leak out from the cells through these transient openings, which were subsequently resealed by replacing the hypotonic solution with an isotonic solution³¹. At the same time, 0.1% w/v sucrose in PBS solution resulted in a sufficiently viscous solution to impede cell rupture or membrane flipping. To maintain the integrity of membrane proteins in n CG, protease inhibitor cocktails were added to CG during the hypotonic solution treatment. As shown in **Figure 3B** and **3C**, more than 70% of the cellular DNA content was removed from 3T3-L1 cells in the process of CG production. This was further confirmed by staining the CG with Hoechst 33342, a fluorescent dye binding to nuclear DNA. The reduced fluorescence reinforces the notion that the DNA was largely removed through the CG production protocol (**Figure 3D**). Noteworthy, the CG remained as spherical vesicles similar to the original cells, which may be indicative that the protocol did not destroy the overall structure of the plasma membranes. The size of CG appeared to be slightly larger than that of the cells, probably due to the interaction between the solutes and phospholipids in the plasma membrane after the physical permeabilization using the hypotonic solution.³²

Removal of cellular and nuclear contents in the CG is an essential step to ensure production reproducibility and safety of the produced mCVTs with regards to minimalizing the amount of genetic material from parent cells that is transferred to target

cells. To confirm that the cellular component of mCVTs is non-living, the non-proliferative nature of CG from 3T3-L1 (CG_{3T3-L1}) was investigated and compared to live cells. The number of CG_{3T3-L1} and 3T3-L1 cells in complete culture medium was plotted against time, as shown in **Figure 3E** for 5 consecutive days. As anticipated, the 3T3-L1 cells showed a typical trend of cell proliferation, indicating normal growth of 3T3-L1 in cell culture. The number of cells started to decline at Day 5 probably due to the fact that the culture became overcrowded. In comparison, the number of CG_{3T3-L1} remained relatively constant throughout the 5 days, suggesting that CG produced through this protocol cannot proliferate and are no longer viable.

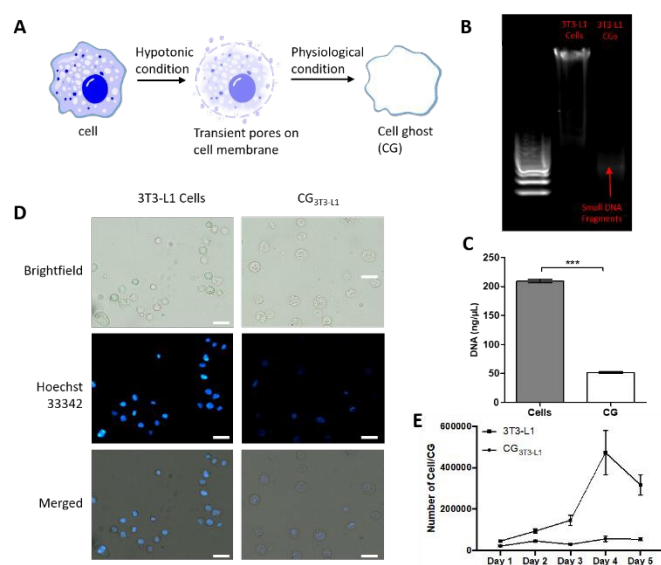


Figure 3. Production and characterization of cell ghosts (CG). **A.** Production scheme of CG by reduction of osmotic pressure. **B.** DNA gel electrophoresis image of 3T3-L1 cells (second lane) and CG (third lane). DNA ladder (first lane) represents 100 bp DNA. **C.** Isolated DNA contents in 3T3-L1 cells and CG, respectively (5×10^6 cells/CG). Data represent means \pm SD ($n=3$), **** $P < 0.0001$. **D.** Fluorescence microscope images of 3T3-L1 cells and CG. Hoechst 33342 was used to stain nuclear contents. Scale bar indicates 100 μ m. **E.** Cell proliferation tests of 3T3-L1 cells and CG. Data represent means \pm SD ($n=3$).

Production and characterization of mCVTs

Similar to liposome production, the lipid thin film method³³ was used to produce mCVTs. DOTAP thin film was hydrated with CG to form large multilamellar liposomes encapsulating CG. This suspension was extruded and plasma membranes from the CG contained DOTAP lipids (as shown in **Figure 4A**). Following a method previously established by our group, it was proven that, upon extrusion, the so-produced mCVTs consisted of the fusion of both lipids and CG membranes (see **Figure 5**). The histogram in **Figure S1** demonstrates the size distribution profile of mCVTs derived from CG_{3T3-L1} (mCVT_{3T3-L1}) and liposomes. Generally, the size distribution of mCVT was similar to liposomes. In order to demonstrate the robustness and versatility of this method across various cell lines, mCVTs were produced using CG derived from different cell lines, such as some HTT cells (i.e. 3T3-L1) and non-cancerous cells (i.e. HEK293 and RAW264.7). As shown in

Figure 4B, the hydrodynamic size and PDI of various mCVTs was typically in the range of 400 nm (with PDI \sim 0.4), similar to DOTAP liposomes (PDI \sim 0.2) produced using the same extrusion method. In **Figure 4C**, similar zeta potential values (around +30mV to +40mV) were observed for different formulations, indicating that this method is generally applicable to CG generated from different cell lines. Interestingly, when mCVT_{3T3-L1} was complexed with pDNA, the size distribution remained similar to that of the non-complexed mCVTs, while the overall zeta potential became more negative (**Figure S2**), which may be indicative for a successful complexation (as the negatively charged pDNA neutralized the positively charged mCVTs).

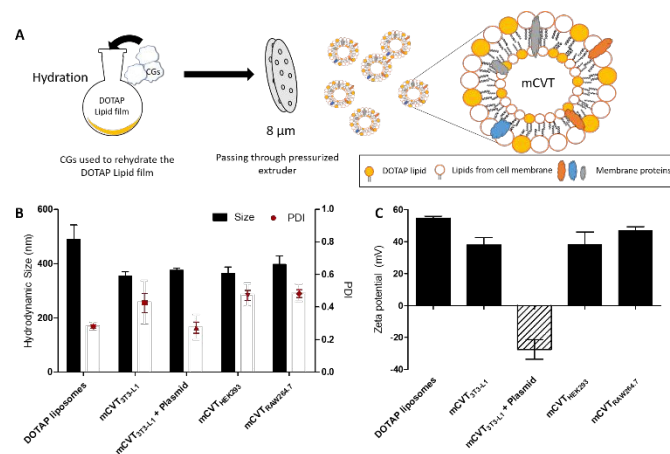


Figure 4. Production and characterization of mCVTs. **A.** Schematic depiction of mCVTs, involving fusion of cell ghosts and liposomes made from thin film hydration. **B.** Hydrodynamic size of different mCVTs and mCVT_{3T3-L1} with plasmid compared to DOTAP liposome. The size of all mCVTs is smaller than 500 nm. **C.** Zeta potential of different mCVTs and mCVT_{3T3-L1} with plasmid compared to DOTAP liposomes. The zeta potential of all mCVTs is above +30 mV, except for mCVT_{3T3-L1} with plasmid DNA (5 μ g), showing negative zeta potential due to the interaction with plasmid DNA. Data represent means \pm SD ($n=3$).

To demonstrate the fusion between the lipids and CG, a Fluorescence Resonance Energy Transfer (FRET) assay was conducted as reported previously.³¹ Two formulations, made from 2.5mg of lipids and either 1×10^7 3T3-L1 CG or 2×10^7 3T3-L1 CG, were used to produce mCVTs and their rhodamine emission intensities were compared against those observed in the DOTAP liposomes. As depicted in **Figure 5A**, in the absence of CG fusion, both the NBD and rhodamine fluorophores in the liposomal formulation were in close proximity, thus DOTAP liposomes alone had the highest observed rhodamine emission intensity at 585 nm. When CG were fused with DOTAP liposomes, cellular membranes from CG were randomly inserted between the two fluorophores, pushing apart the NBD and rhodamine fluorophores, thus resulting in a reduction of the observed emission intensity at 585 nm. Consistent with this observation, the emission intensity of rhodamine decreased further when twice the amount of CG was used for mCVTs production (**Figure 5A**). Altogether, the results confirmed the fusion between cationic lipids and CG, rather than the encapsulation of CG into liposomes or *vice versa*.

mCVTs produced from various amount of CG (**Table S3**) were characterized by measuring the hydrodynamic size, zeta potential and protein concentration. Increasing the CG quantity in the mCVTs formulation reduced the hydrodynamic size of the prepared mCVTs (from half amount, i.e. mCVT 0.5x, to full amount, i.e. mCVT 1x), till a hydrodynamic mean size of about 400 nm was observed (in mCVT 1x). When more CG were added (double amount, i.e. mCVT 2x), the hydrodynamic size of the mCVTs increased, possibly due to excess of CG leading to a higher propensity to aggregate (**Figure 5B**). Moreover, the size of mCVTs was unaffected after the complexation, suggesting that the plasmid DNA was complexed to mCVTs by electrostatic interactions and did not result in aggregation of mCVTs. The zeta potential was generally positive due to the incorporation of cationic DOTAP. Since CG have a negative zeta potential as inherited from the parent cells, the overall positive charge of mCVTs confirms the successful incorporation of DOTAP into the cell membrane of CG. The positive zeta potential decreased when a higher concentration of CG_{3T3-L1} was used (e.g. with mCVT 2x), providing additional evidence that DOTAP lipids were indeed inserted into the plasma membranes. The greatest decrease in zeta potential was observed when 1×10^7 3T3-L1 cells were used to interact with DOTAP lipids (mCVT 1x (cells)), showing a negative zeta potential. These decreases in zeta potential are likely due to the increasing amounts of proteinaceous materials or other cellular components (e.g. nuclear content) from the cells relative to the amounts of cationic lipids (**Figure 5C**). Once again, this highlights the importance of producing CG prior to mCVTs production as the interaction of cationic lipids with cells can immediately result in formation of aggregates instead of mCVTs (**Figure 5D**). CG were also not a suitable construct for transfection as no transfection was observed with cells treated with the CG alone. This is likely due to the inability of CG (which have a negative zeta potential) to complex with plasmid (**Figure 5E**). Interestingly, mCVT 2x did not show a higher protein concentration or better transfection performance compared to mCVT 1x, despite the fact that only half of the amount of CG was used for the production of the latter preparation (**Figure 5D and 5E**). One plausible explanation is a lower protein yield in case of mCVT 2x due to the formation of protein aggregates during the extrusion process, as the higher amount of CG leads to much more clogging of the extrusion membranes, which translated into substantially more protein loss.

Cell Ghosts (CG) are required for production of mCVTs

To investigate the effect of cellular contents on the production yield and the transfection efficiency of a model gene (enhanced green fluorescence protein, eGFP), 1×10^7 of 3T3-L1 cells were used for mCVTs production and compared with mCVTs produced from CG. mCVT 1x (cells) showed sub-micron size with a negative net zeta potential and high protein concentration, which did not contribute to the *in vitro* gene delivery to cells. Similar to mCVT 2x from CG, protein aggregation happened during the production process with mCVT 1x (cells), and the bulk of proteins from both plasma membranes and intracellular contents hindered the interaction

with cationic DOTAP, resulting in a net negative charge with no transfection. Conversely, the mCVT 1x (from CG rather than whole cells) showed to be a suitable combination of cellular membrane (and proteins) and lipid charges to form mCVTs, which translated into a high transfection efficiency and viability for target HTT cells (**Figure 5E, 5F**) in comparison to DOTAP liposomes or CG alone. Altogether, these results emphasize the importance of the cellular component from the CG in providing plasmid transfection in HTT cells, as well as in masking the cytotoxicity of cationic liposomes.

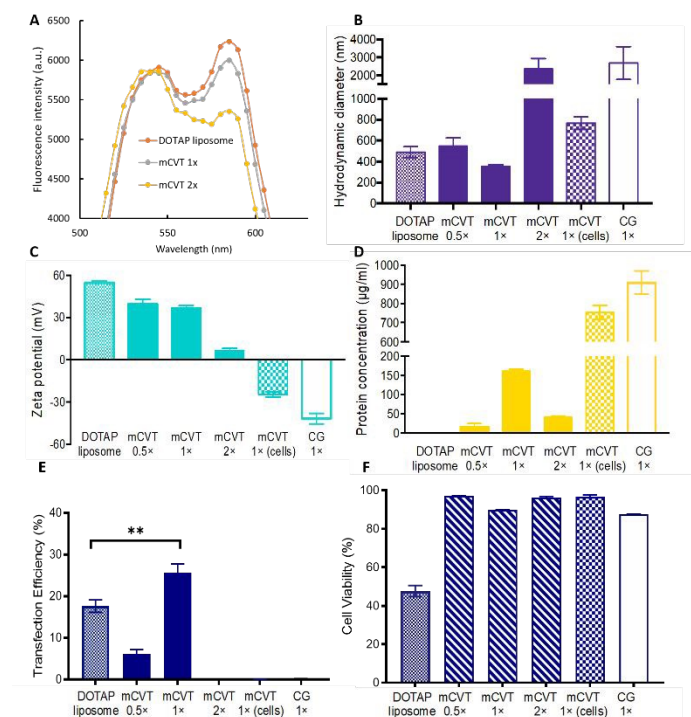


Figure 5. Proof-of-fusion of mCVT_{3T3-L1}. Fluorophores are in close proximity in liposomes, while fusion of DOTAP liposomes with CG increases the distance between two fluorophores. **A.** FRET analysis of NBD and rhodamine fluorophores of DOTAP liposome and mCVTs. Upon exciting at 460 nm, the energy of emission spectrum for NBD excited rhodamine, resulting in the emission of rhodamine at 585 nm. The emission of rhodamine was higher compared to mCVT 1x and mCVT 2x, suggesting fusion between DOTAP liposomes and CG_{3T3-L1}. **B.** Hydrodynamic size **(B)** Zeta potential **(C)** and protein concentrations **(D)** of DOTAP liposome, mCVT 0.5x, mCVT 1x, mCVT 2x, DOTAP + cells and CG. **E.** Transfection efficiency **(E)** and cell viability **(F)** of 3T3-L1 cells upon transfection. n=3, **p<0.01.

Presence of plasma membrane components improves the cellular uptake of mCVTs by HTT cells

In order to understand how plasma membrane components in mCVTs affect uptake and transfection of target HTT cells, cyanine 5.5 (Cy5.5) was used to label the proteins in the outer cellular membrane of CG_{3T3-L1}, which were used to produce mCVT_{3T3-L1}; and NBD labelled DOTAP lipids were used for the formation of the cationic liposomes.

Figure 6A shows the fluorescence profiles by flow cytometry of CG_{3T3-L1}, DOTAP liposomes and mCVT_{3T3-L1} after labelling. The presence of both NBD and Cy5.5 in mCVT_{3T3-L1} indicates the incorporation of both DOTAP lipids and 3T3-L1 cell membrane components in mCVT_{3T3-L1}. This is in agreement with our FRET

results (**Figure 5A**), confirming that mCVT_{3T3-L1} is a hybrid system consisting of these two components.

Subsequently, the fluorescence-labelled DOTAP liposomes, CG_{3T3-L1} and mCVT_{3T3-L1} were added to 3T3-L1 cells to study their uptake profile (**Figure 6B**). DOTAP liposomes were taken up by the cells at a slower rate as compared to CG_{3T3-L1}, as after 1 hour of incubation with liposomes, only a small percentage of cells were positive for NBD (< 10%), whereas > 90% of cells were positive for Cy5.5 after 1 hour of incubation with CG_{3T3-L1}. For the cellular uptake of mCVT_{3T3-L1}, it appeared that our mCVTs (hybrid system) were being taken up by cells at a rate comparable to CG_{3T3-L1} (plasma membrane particles), and faster than DOTAP liposomes (fully synthetic lipid vesicles). Thus, as mCVTs represent a hybrid system of cationic lipids and cell membrane components, we postulate that they facilitate the internalization of plasmid and transfection through the following mechanisms: (1) complexation of plasmid DNA with cationic lipids allows mCVTs to carry the plasmid DNA into 3T3-L1 cells, as naked plasmid is unable to enter cells on its own; (2) the positive charges of cationic lipids facilitate the interaction of mCVTs with the negatively charged cell membranes of 3T3-L1 cells; (3) the surface properties (from plasma membrane components) on mCVTs further promote the interaction with target cells; and (4) efficient cellular uptake of membrane components of mCVTs (as indicated by the efficient uptake of CG) facilitates the delivery of plasmid DNA into the cells. Taken together, the presence of membrane components in mCVTs facilitates transfection of plasmid DNA into HTT cells.

mCVTs display superior transfection capability and reduced cytotoxicity

Having characterized mCVT_{3T3-L1} and demonstrated their cellular uptake, we next investigated plasmid delivery efficiency of mCVT_{3T3-L1} using 3T3-L1 cells and eGFP as model plasmid. Transfection efficiency of mCVT_{3T3-L1} was compared with a widely used commercial transfection reagent, Lipofectamine 3000 (LF3000). When 5 µg of the total plasmids were used, the maximal transfection efficiency (expressed as number of cells transfected) of 3T3-L1 cells, as a HTT cell line, using LF3000 was about 15% at 10 µl (**Figure 7A**). Higher amounts of LF3000 resulted in a lower transfection efficiency attributable to a higher cytotoxicity (**Figure 7B**). In comparison, the maximal transfection efficiency of 3T3-L1 cells by mCVT_{3T3-L1} (equivalent to 1.25 mg/ml of DOTAP lipids, and about 80 µg/ml of CG proteins) reached 30% at 80 µl. Because mCVT_{3T3-L1} started to induce detectable cytotoxicity in 3T3-L1 cells at the amount of 80 µl, we considered 40 µl of mCVT_{3T3-L1} as the optimal transfection dose that reached significantly higher transfection efficiency (25±3.1% for mCVT_{3T3-L1} vs. 15±0.3% for LF3000, *p*<0.001) with marginal cytotoxicity.

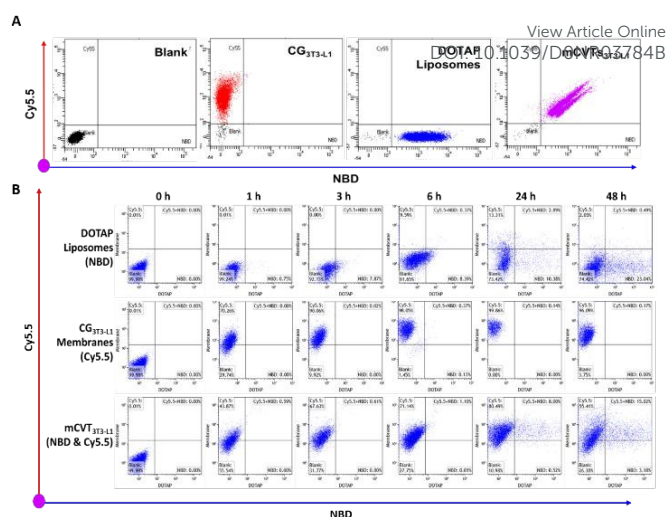


Figure 6. Uptake of mCVT_{3T3-L1}, DOTAP liposome and CG_{3T3-L1} by 3T3-L1 cells. Fluorescently labelled components (NBD-DOTAP and Cy5.5-CG_{3T3-L1}) were used for production of mCVTs. **A.** Dot plots of CG_{3T3-L1}, DOTAP liposome and mCVT_{3T3-L1} for flow cytometry analysis. Black: Blank; Red: Cy5.5; Blue: NBD; Purple: both Cy5.5 and NBD. **B.** Uptake of DOTAP liposome, CG_{3T3-L1} and mCVT_{3T3-L1} into 3T3-L1 cells. Analysis was done for 1, 3, 6, 24 and 48 h. Blank: Green; Red: Cy5.5; Blue: NBD; Purple: both Cy5.5 and NBD.

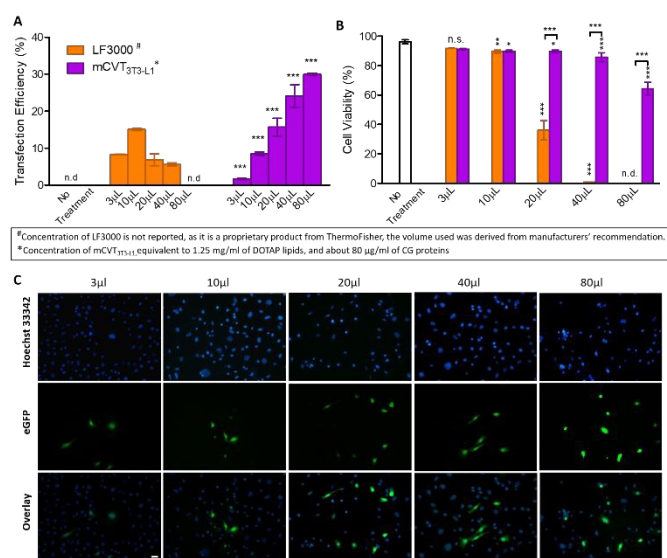


Figure 7. Transfection efficiency and cytotoxicity of mCVTs. **A.** Transfection efficiency of LF3000 and mCVT_{3T3-L1} into 3T3-L1 cells. 5 different amounts (3 µl, 10 µl, 20 µl, 40 µl and 80 µl) of LF3000 were used to transfect 5 × 10⁵ number of 3T3-L1 cells, with the same corresponding doses of mCVTs chosen for the same transfection experiments for comparison. Cells without treatment were used as the control. *** *P* < 0.001 (for comparison to respective amount of LF3000) **B.** Cell viability of 3T3-L1 after transfection using LF3000 and mCVT_{3T3-L1}. Cells without treatment were used as the control. n.d. and n.s. indicated no detectable transfection and not significant respectively. Data represented means ± SD (*n*=3), * *P* < 0.05 ** *P* < 0.01, *** *P* < 0.001. **C.** Images of transfection cells using different doses of mCVT_{3T3-L1}. The transfected cells expressed eGFP (green). Hoechst 33342 (blue) was used to visualize the nucleus. Scale bar represents 100 µm.

The mCVT_{3T3-L1} transfected cells were also monitored using fluorescence microscopy imaging. The cell nuclei were visualized using Hoechst 33342 (**Figure 7C**). In agreement with **Figure 7A**, an increase in intracellular eGFP signal was observed when higher doses of mCVT_{3T3-L1} were used. Overall, mCVT_{3T3-L1} were about twice as effective as LF3000 with regards to transfection efficiency and yet exhibited only moderate cytotoxicity compared to LF3000, indicating their significant advantage for the transfection of HTT cells.

Next, we explored whether the use of mCVTs for transfection was associated with some extent of specificity towards the target cells. mCVTs were produced using CG from different cell types, including 3T3-L1 cells, HEK293 cells and RAW264.7 macrophages, and were tested for their ability to transfect our model gene into 3T3-L1 cells. LF3000, CG and DOTAP liposomes were used as controls (**Figure 8A**). The dosage of LF3000 (10 μ l), DOTAP liposomes (40 μ l) and mCVTs (40 μ l) were selected based on maximal transfection efficiency and minimal cytotoxicity results (see **Figure 8B**).

Cells treated with CG_{3T3-L1} and plasmid showed no detectable transfection due to the inability of anionic CG to complex with plasmid DNA. Although it was shown that CG could be taken up by cells without complexation, CG were not able to facilitate the entry of plasmid DNA into 3T3-L1 cells. DOTAP liposomes exhibited a similar transfection efficiency as LF3000 at a dose of 40 μ l but they caused high cytotoxicity due to their highly cationic nature. Conversely, all mCVTs displayed minimal cytotoxicity towards 3T3-L1 cells and they showed either comparable (in the case of mCVT_{HEK293} and mCVT_{RAW264.7}) or superior transfection efficiency (in the case of mCVT_{3T3-L1}) when compared to the commercial LF3000 and DOTAP liposomes. This is likely due to the presence of plasma membrane components (from CG) that rendered mCVTs less cytotoxic. This again highlights the advantage of using mCVTs as transfection agent.

Among the three different mCVTs, mCVT_{3T3-L1}, mCVT_{HEK293} and mCVT_{RAW264.7}, mCVT_{3T3-L1} showed the highest transfection efficiency (**Figure 8A**). Since not all mCVTs (dose dependent transfection efficiency of different mCVTs and respective cell viability can be found in **Figure S3**) were able to transfect 3T3-L1 to a similar extent, we reasoned that the presence of plasma membrane components of 3T3-L1 cells may have induced specificity towards 3T3-L1 cells. This homotypic targeting behaviour suggests that, when formulating mCVTs, the selection of the appropriate cell line is crucial in improving the overall transfection efficiency.

In this study, 3T3-L1 cells, a murine preadipocyte cell line, were used as a quintessential example of HTT cells. mCVTs produced from 3T3-L1 CG were able to improve the transfection efficiency to approximately 30% while maintaining nearly 90% of cell viability upon transfection. This is in contrast with other commercial transfection reagents (for example LF3000), which have limited transfection efficiency (\approx 15%), and high cytotoxicity when used at similar concentrations. Furthermore, mCVTs were found to maintain a relatively high transfection efficiency (\approx 20%), even under serum condition (**Figure S4**). This highlights the additional advantage of using mCVTs as potential *ex vivo* or *in vitro* transfection reagent, as serum can be added during the process to nourish and protect the sensitive cell lines (i.e. primary cells).

By studying mCVTs and the extent of cellular uptake when different parent cells were used, we identified several key insights. First, the transfection efficiency of non-viral-based gene delivery systems is typically influenced by the zeta potential of the overall complexes, comprising of vector and genetic load. Data from the literature recommend that a higher positive zeta potential is preferred to increase gene transfection,^{34,35} but that cell viability is strongly limiting the degree of transfection. Surprisingly, we did not observe such a correlation for mCVTs. In our study, mCVTs had an overall negative zeta potential between -30 to -40 mV after complexing with the plasmid, and yet they were able to successfully transfect the HTT cell lines. We speculate that the plasma membrane components present in mCVTs likely contribute to improving the transfection efficiency in case of HTT cells. Second, the use of different precursor cells in the production of mCVTs affects the overall transfection efficiency and cytotoxicity profiles. Interestingly, we observed a certain degree of “like-to-like” selectivity of mCVTs_{3T3-L1} towards 3T3-L1 cells, as mCVTs_{3T3-L1} transfect 3T3-L1 cells more efficiently than the mCVTs produced from other cell lines (similar results were obtained when HEK cells were transfected with mCVT_{HEK293} (**Figure S5**)). We postulate that mCVTs are preferentially taken up by the cells of the same origin to some degree, likely due to the conservation of surface recognition moieties during the production of mCVTs. These fingerprints are later used as recognition moieties for homotypic targeting. A similar phenomenon was reported for other cell-based DDS (i.e. exosomes). Sancho-Albero *et al.*³⁶ demonstrated a selective uptake of exosomes derived from mesenchymal stem cells

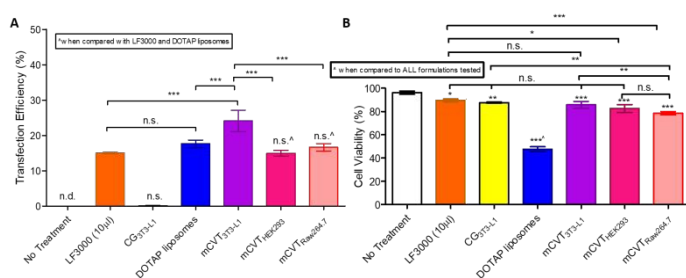


Figure 8. Transfection efficiency and cytotoxicity of different mCVTs on 3T3-L1 cells after 24 hours incubation. **A.** Transfection efficiency of LF3000 (10 μ l), CG_{3T3-L1}, DOTAP liposomes (40 μ l) and different mCVTs (40 μ l) on 3T3-L1 cells. **B.** Cell viability of 3T3-L1 cells after transfection with LF3000, CG_{3T3-L1}, DOTAP liposomes and different mCVTs. n.d. and n.s. indicated no detectable transfection, and not significant respectively. Data represented means \pm SD ($n=3$), * $P < 0.05$, ** $P < 0.01$ and *** $P < 0.001$.

(MSC) by other MSCs in a co-culture system. Likewise, the results reported by Hazan-Halevy *et al.*³⁷ indicated that exosomes derived from B lymphocytes were efficiently internalized by other B lymphocytes and not by other immune cells. All these studies suggest a role for surface recognition moieties in specifically recognising the cell types used for the production of the mCVTs.

However, it is noteworthy that the nature of the selective targeting of cell-based nanovesicles (e.g. exosomes) is still largely debatable: some studies (like those reported above) have demonstrated homotypic targeting, while others highlighted that cell-based vesicles could be used to target different target cells, such as nanovesicles derived from U937 monocytes showing targeting properties towards inflammatory or cancer tissue.^{25,26}

Therefore, the selective nature of mCVTs may not be solely towards precursor cells. Instead, we found that it was possible to use mCVT_{3T3-L1} to transfect other HTT cells such as keratocytes (HaCaT) (Figure S6). One plausible explanation is that mCVTs possess receptor/ligand molecules, acquired from precursor cells, which are responsible for cell-cell interactions. In this regard, 3T3-L1 cells are known to express a basal level of inflammatory markers such as MCP-1 and ICAM-1, and other integrin molecules which are responsible for cell adhesions.³⁸ When passing down to mCVT_{3T3-L1}, these cell surface markers may hence enable mCVT_{3T3-L1} to interact with other cells with corresponding receptors/ligands, leading to better uptake of mCVTs and, in this case, improved transfection. Our results thus lay the groundwork to further explore the possible interactions between membrane proteins across different cell types.

Overall, we have demonstrated here the advantages of the presence of plasma membrane components in mCVTs in order to enhance selectivity and cellular uptake, as well as to increase the overall transfection efficiency of HTT cells. Furthermore, we demonstrated the “tunability” of producing mCVTs from various precursor cells according to the desired transfection outcome. As a hybrid system of plasma membrane components and lipids, various lipids can be added to or replace the cationic lipids, in order to confer new functionalities to mCVTs. This flexibility could be further exploited for the addition of fluorophores conjugated to lipids for visualization or tracking purposes, or for incorporation of magnetic nanoparticles for cell manipulation or even for attaching additional cell targeting ligands to achieve improved target cell selectivity.

Conclusions

In this study, we have established a protocol for CG production by varying the osmotic pressure, where CG were shown to be non-proliferative, with most nuclear material removed. Using cationic lipids and CG, we constructed a hybrid plasmid delivery system comprising both plasma membrane and cationic lipid components, which we termed mCVTs. We have also shown that mCVTs have similar characteristics and ability to complex with plasmid DNA as cationic liposomal systems. Furthermore,

mCVTs are superior to cationic liposomes and the commercial transfection reagent LF3000, in terms of transfection efficiency and cytotoxicity. Due to the chimeric nature of mCVTs, comprising both plasma membrane components and cationic lipids, mCVTs were more efficiently taken up by target cells. Furthermore, mCVTs, being cell-derived formulations, confer flexibility in changing the cell-based components depending on the target cells or tissues and do not require heavy capital investment into sophisticated pieces of equipment. Hence, the whole process of production of mCVTs is ideal for quick in-house testing and also amenable for up-scaling. Taken together, the successful development of mCVTs provides a basis for the *in vitro* delivery of plasmids and other nucleic acids, especially for HTT cells, paving the way for gene delivery into a broader range of cell types.

Material & Methods

Materials

1,2-dioleoyl-3-trimethylammonium-propane (DOTAP), 1-oleoyl-2-[6-[(7-nitro-2-1,3-benzoxadiazol-4-yl)amino]hexanoyl]-3-trimethylammonium propane (NBD-DOTAP) and 1,2-dioleoyl-*sn*-glycero-3-phosphoethanolamine-N-(lissamine rhodamine B sulfonyl) (Liss Rhodamine) were purchased from Avanti Polar Lipids (Alabama, United States) and kept at -20 °C. Chloroform was used to dissolve the lipids. Lipofectamine 3000 (LF3000), Hoechst 33342 nuclear stain, trypan blue stain solution (Hyclone™) and Bicinchoninic Acid (BCA) protein assay and eGFP expressing plasmid was purchased from ThermoFisher Scientific (Massachusetts, United States). 8 µm membrane filters were purchased from Merck Millipore (Massachusetts, United States) and used as supplied. Cyanine 5.5 *N*-hydroxysuccinimide (NHS) monoester dyes were purchased from Kerfast (Boston, United States) and used as per manufacturer's recommendations. Polybead® Microspheres 10 µm were obtained from Polysciences, Inc. (Pennsylvania, United States).

Cell Culture

3T3-L1 cells, a fibroblast-like mouse preadipocyte cell, were cultured in Dulbecco's High Glucose Modified Eagles Medium/high glucose (DMEM/high glucose, HyClone™) containing 10% v/v foetal bovine serum (FBS, HyClone™) and 1% v/v MEM non-essential amino acid solution (NEAA, Life Technologies, California, United States). HEK293 cells (human embryonic kidney cells) and RAW264.7 cells (mouse macrophages), were cultured in DMEM/high glucose supplemented with 10% FBS. Cells were cultured in a humidified incubator maintained at 5% CO₂ and 37 °C.

Production of Cell Ghosts (CG)

CG were produced *via* changing osmotic pressure using neutral phosphate buffer saline (PBS, pH 7.4) with 0.1% w/v sucrose (PBS/sucrose). Cells were harvested and incubated in 0.25x PBS (diluted by sterile MilliQ water) at room temperature for 24 hours, followed by incubating CG in 0.1x PBS/sucrose. CG were harvested in 1x PBS/sucrose and stored at 4°C for further experiments. Protease inhibitor cocktail (Abcam, Cambridge,

United Kingdom) was added throughout the experiment. The presence of nuclei material was checked by staining CG with Hoechst 33342. And remaining DNA material in CG was isolated using Qiagen DNeasy Blood and Tissue Kit (Hilden, Germany) for DNA gel electrophoresis according to manufacturer's instruction.

Cell ghosts (CG) viability test

The same number of 3T3-L1 cells and CG were prepared in DMEM high glucose culture medium with 10% FBS and 1% NEAA. The samples were collected, and total cell or CG count were determined by staining with trypan blue solution and counting *via* a cell counting chamber on days 1, 2, 3, 4 and 5.

Production of mCVTs

2.5 mg of DOTAP was dissolved in 1 ml of chloroform in a 5 ml round bottom flask. A dry thin lipid film was formed on the round bottom flask under vacuum at 35°C for at least 1 hour. CG suspended in PBS were added into the dry lipid film and were subsequently agitated by brief sonication. The mixture was passed through 8 µm membrane filters using a liposome extruder (Genizer™, California, United States) for 3 times. Samples were taken for the measurement of hydrodynamic size, Poly Dispersity Index (PDI) and zeta potential. Liposomes used in this study were produced, in the same manner, with same formulation without adding CG.

mCVTs characterization

All hydrodynamic size and zeta potential measurements were determined at 25°C by Dynamic Light Scattering (DLS) on a Zetasizer Nano ZS (Malvern, UK). Characterization data were analyzed using software supplied with the Zetasizer and the results from at least three different measurements per run for more than 10 runs carried per sample were presented. Protein concentration was determined using (bicinchoninic acid assay) BCA protein assay for each sample.

Proof-of-fusion assay

A Fluorescence Resonance Energy Transfer (FRET) assay was conducted. Liposome and mCVT_{3T3-L1} were prepared as described previously³¹ with the addition of 0.1 mol% NBD-DOTAP and lissamine rhodamine lipid. To prove the incorporation of cell membrane components from CG_{3T3-L1} to mCVT_{3T3-L1}, various concentrations of CG (0.5 × 10⁷, 1 × 10⁷ and 2 × 10⁷ of CG_{3T3-L1}) were used to produce mCVT_{3T3-L1} while the amount of DOTAP used remain unchanged. The liposome and mCVTs were then measured at 460 nm using a microplate reader and the emission spectra from 300 nm to 700 nm were recorded.

Transfection of 3T3-L1 using LF3000 and mCVTs

3T3-L1 cells with 60-70% confluency were ready for transfection. 3T3-L1 cells were transfected with plasmid DNA containing pEZX-MR04 (miRNA overexpression scrambled control clone, GeneCopoeia) using Lipofectamine 3000 according to the manufacturer's instruction. For transfection of 3T3-L1 using mCVTs, plasmid DNA (fixed complexation ratio of

5 µg/40 µl) was incubated with various concentration of mCVTs at room temperature for 10 minutes, followed by incubating the cells with serum-free medium. The final amount of pDNA used in corresponding experiments (LF3000 and mCVTs) were the same. After incubating in serum free medium for 6 hours, the cells were rinsed with PBS and cultured under complete medium with 1% Active Ingredient-Active Ingredient (Hyclone, Utah, United States). Propidium iodide (PI) was added to determine cell viability, the viable cells with intact cell membrane will be able to exclude the dye, whereas, dead/dying cells will be stained. The percentage transfection efficiency and cell viability were determined by detecting eGFP (reporter of the plasmid) and PI (as cell viability dye), respectively, using BD LSR Fortessa flow cytometry analyser (New Jersey, United States). The untreated cells were used as control. All fluorescence images were taken using Olympus BX51 fluorescence microscope (Tokyo, Japan). O

Flow cytometry analysis of mCVTs and liposomes

Fluorophore labelled DOTAP liposomes were made by incorporating 0.1 mol% NBD-DOTAP in the lipid formulation of DOTAP liposome. Cyanine 5.5 (Cy5.5) *N*-hydroxysuccinimide (NHS) monoester was used for the labelling of membrane protein of CG. The fluorophore labelled CG were then used to produce mCVTs. DOTAP liposomes and mCVTs were adsorbed on to the 10 µm polybeads (through electrostatic interaction) by incubating these vesicles with the beads at 4°C overnight under gentle agitation. The fluorophore intensity was analyzed by using BD LSR Fortessa flow cytometry analyser.

Cellular Uptake of mCVTs

Fluorophore labelled DOTAP liposomes, CG and mCVTs were produced as described above for the cellular uptake experiment. 3T3-L1 cells with 60-70% confluency were incubated with fluorophore labelled DOTAP liposomes, CG and mCVTs. The fluorophore intensity was analyzed after 1, 3, 6, 24 and 48 hours of incubation by using BD LSR Fortessa flow cytometry analyser.

Statistical Analysis

All statistical analyses were performed using GraphPad Prism 7 using one-way ANOVA, followed by Bonferroni *post hoc* test. *P* < 0.05 was deemed significant.

Conflicts of interest

The authors declare the following competing financial interest(s): Wei Jiang Goh, Shui Zou and Giorgia Pastorin are the inventors of a patent relating to the production and use of CVTs as drug delivery vehicles (PCT/SG2018/050376; Title: Biomolecular Composites Comprising Modified Cell Ghosts).

Acknowledgements

This work was supported by National University of Singapore (NUS), Department of Pharmacy and by the funding provided by MOE (grant number R-148-000-296-114) and partly supported by the Singapore NMRC (NMRC/OFYIRG/0081/2018 to J.W.W., G.S., B.C. and G.P.). The authors would like to acknowledge A/P Gigi

Chiu and A/P Ang Wee Han in NUS for the use of cells and equipment.

Reference

- 1 S. Yamano, J. Dai and A. M. Moursi, *Mol. Biotechnol.*, 2010, **46**, 287–300.
- 2 T. K. Kim and J. H. Eberwine, *Anal. Bioanal. Chem.*, 2010, **397**, 3173–3178.
- 3 J. F. Tan, T. A. Hatton, K. C. Tam and H. P. Too, *Biomacromolecules*, 2007, **8**, 448–454.
- 4 A. K. Pannier and L. D. Shea, *Mol. Ther.*, 2004, **10**, 19–26.
- 5 T. Niidome and L. Huang, *Gene Ther.*, 2002, **9**, 1647–1652.
- 6 Y. Seow and M. J. Wood, *Mol. Ther.*, 2009, **17**, 767–777.
- 7 G. D. Schmidt-Wolf and I. G. H. Schmidt-Wolf, *Trends Mol. Med.*, 2003, **9**, 67–72.
- 8 E. Jubeli, W. P. D. Goldring and M. D. Pungente, in *Methods in Molecular Biology*, Humana Press Inc., 2016, **1445**, 19–32.
- 9 P. Midoux, G. Breuzard, J. Gomez and C. Pichon, *Curr. Gene Ther.*, 2008, **8**, 335–352.
- 10 F. Prüfert, S. Bonengel, S. Köllner, J. Griesser, M. D. Wilcox, P. I. Chater, J. P. Pearson and A. Bernkop-Schnürch, *Nanomedicine*, 2017, **12**, 2713–2724.
- 11 D. Lou and W. M. Saltzman, *Nat. Biotechnol.*, 2000, **18**, 33–37.
- 12 S. L. Holmen, M. W. Vanbrocklin, R. R. Eversole, S. R. Stapleton and L. C. Ginsberg, *Vitr. Cell. Dev. Biol. - Anim. J. Soc. Vitr. Biol.*, 1995, **31**, 347–351.
- 13 P. Washbourne and A. K. McAllister, *Curr. Opin. Neurobiol.*, 2002, **12**, 566–73.
- 14 X. Han, Z. Liu, M. C. Jo, K. Zhang, Y. Li, Z. Zeng, N. Li, Y. Zu and L. Qin, *Sci. Adv.*, 2015, **1**, e1500454.
- 15 G. Kilroy, D. H. Burk and Z. E. Floyd, *PLoS One*, 2009, **4**, e6940.
- 16 S. M. Tietz and M. Berghoff, *J. Biomol. Tech.*, 2012, **23**, 47–50.
- 17 N. Zhao, J. Qi, Z. Zeng, P. Parekh, C. C. Chang, C. H. Tung and Y. Zu, *J. Control. Release*, 2012, **159**, 104–110.
- 18 J. Yin, Z. Ma, N. Selliah, D. K. Shivers, R. Q. Cron and T. H. Finkel, *J. Immunol. Methods*, 2006, **312**, 1–11.
- 19 L. Novo, E. Mastrobattista, C. F. VanNostrum, T. Lammers and W. E. Hennink, *Expert Opin. Drug Deliv.*, 2015, **12**, 507–512.
- 20 F. J. Ruiz-Ojeda, J. Plaza-Díaz, M. J. Sáez-Lara and A. Gil, *Adv. Nutr.*, 2019, **10**, S31–S48.
- 21 O. Gresch and L. Altrogge, *Methods Mol. Biol.*, 2012, **801**, 65–74.
- 22 A. Zotova, E. Lopatukhina, A. Filatov, M. Khaitov and D. Mazurov, *Viruses*, 2017, **9**, 325.
- 23 H. F. Teixeira, F. Bruxel, M. Fraga, R. S. Schuh, G. K. Zorzi, U. Matte and E. Fattal, *Int. J. Pharm.*, 2017, **534**, 356–367.
- 24 W. J. Goh, C. K. Lee, S. Zou, E. C. Woon, B. Czarny and G. Pastorin, *Int. J. Nanomedicine*, 2017, **12**, 2759–2767.
- 25 S. C. Jang, O. Y. Kim, C. M. Yoon, D.-S. Choi, T.-Y. Roh, J. Park, J. Nilsson, J. Lötvall, Y.-K. Kim and Y. S. Gho, *ACS Nano*, 2013, **7**, 7698–7710.
- 26 W. J. Goh, S. Zou, W. Y. Ong, F. Torta, A. F. Alexandra, B. M. Schiffelers, G. Storm, J.-W. Wang, B. Czarny and G. Pastorin, *Sci. Rep.*, 2017, **7**, 14322. DOI: 10.1038/s41598-017-07918-8
- 27 L. Kaneti, T. Bronshtein, N. Malkah Dayan, I. Kovregina, N. Letko Khait, Y. Lupu-Haber, M. Fliman, B. W. Schoen, G. Kaneti and M. Machluf, *Nano Lett.*, 2016, **16**, 1574–1582.
- 28 H. M. Byun, D. Suh, H. Yoon, J. M. Kim, H. G. Choi, W. K. Kim, J. J. Ko and Y. K. Oh, *Gene Ther.*, 2004, **11**, 492–496.
- 29 T. Li, H. Dong and R. Mo, *Nano Res.*, 2018, **11**, 5240–5257.
- 30 E. V. Batrakova and M. S. Kim, *J. Control. Release*, 2015, **219**, 396–405.
- 31 W. J. Goh, S. Zou, B. Czarny and G. Pastorin, *Nanoscale*, 2018, **10**, 6812–6819.
- 32 W. F. Wolkers, Harrië Tte Oldenhof, F. Tang, J. Han, J. Bigalk and H. Sieme, *Langmuir*, 2019, **35**, 7520–7528.
- 33 C. Has and P. Sunthar, *J. Liposome Res.*, 2019, 1–30.
- 34 L. Huang and S. Li, *Nat. Biotechnol.*, 1997, **15**, 620–621.
- 35 M. Agirre, J. Zarate, G. Puras, E. Ojeda and J. Luis Pedraz, *Drug Deliv.*, 2015, **22**, 100–110.
- 36 M. Sancho-Albero, N. Navascués, G. Mendoza, V. Sebastián, M. Arruebo, P. Martín-Duque and J. Santamaría, *J. Nanobiotechnology*, 2019, **17**, 16.
- 37 I. Hazan-Halevy, D. Rosenblum, S. Weinstein, O. Bairey, P. Raanani and D. Peer, *Cancer Lett.*, 2015, **364**, 59–69.
- 38 A. L. Dordevic, N. Konstantopoulos and D. Cameron-Smith, *PLoS One*, 2014, **9**, e99382.
- 39 D. E. Murphy, O. G. deJong, M. Brouwer, M. J. Wood, G. Lavieu, R. M. Schiffelers and P. Vader, *Exp. Mol. Med.*, 2019, **51**, 1–12.

A new hybrid gene delivery platform (mCVT) remarkably improves the transfection efficiency of Hard-To-Transfect (HTT) cells, while preserving excellent biocompatibility.

View Article Online
DOI: 10.1039/D0NR03784B

